


Flexible Tip Bougie™

Articulating Tip, Easy Slide Control

- Increases chance of successful intubation
- Articulating tip flexes into best position to pass through vocal cords
- Easy one-handed operation
- High visibility tip glows in UV light

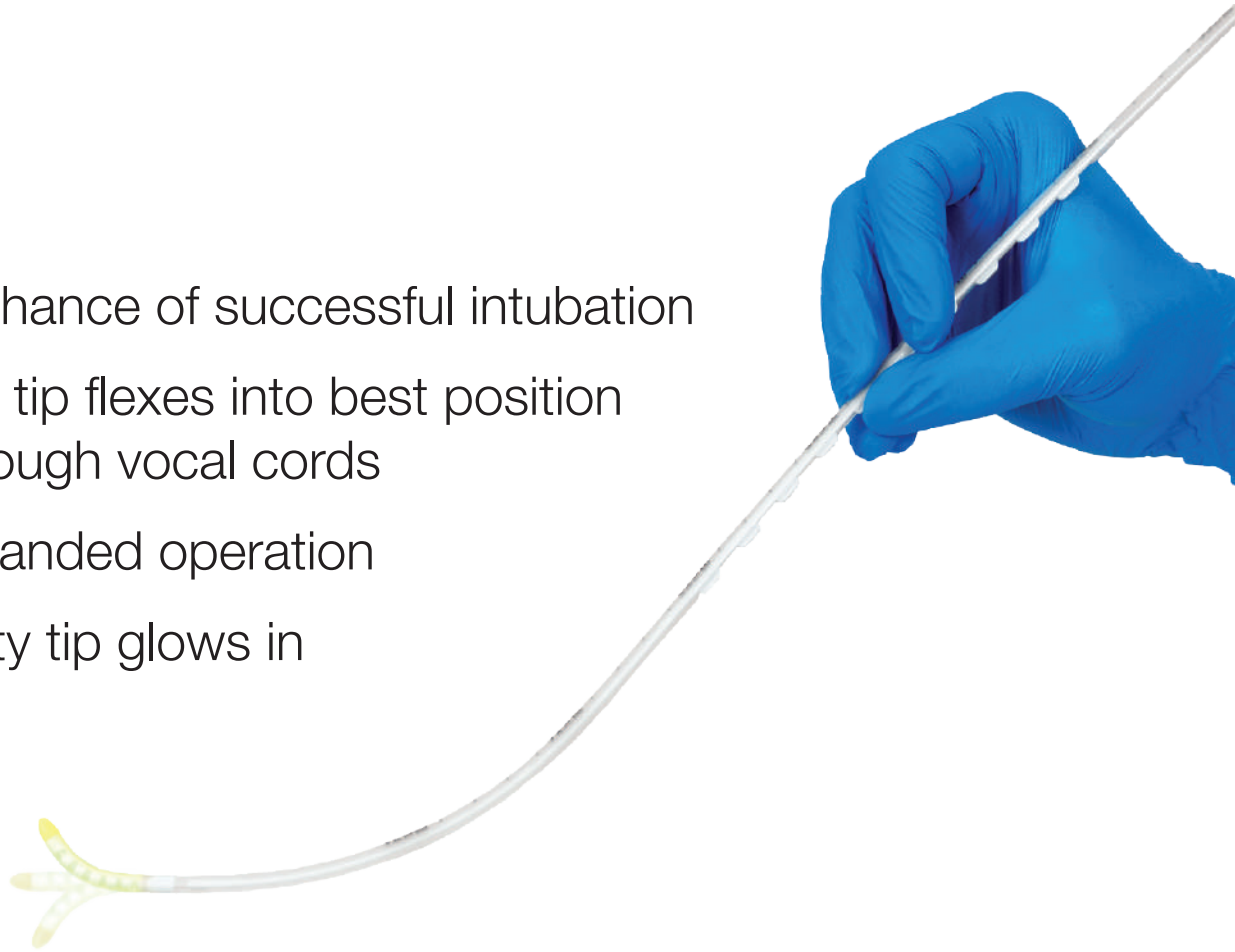

Complete Control at Your Fingertips

The Flexible Tip Bougie™ is a universal, disposable, and articulating bougie for use with adult oral endotracheal tubes with an inner diameter of 7mm and higher. Its one-handed operation makes it an ideal aid when establishing an airway. The tip of the device can easily be maneuvered up and down by manipulating the device's slider tabs while inserting the fluorescent tip through the vocal cords.

Flexible Tip Bougie™

P/N: BC-01

Product Highlights

50cm depth

Fluorescent tip

Latex-free



This universal bougie is ideal for any caregiver's bag

- Controllable tip to facilitate placement through the vocal chords
- Centimeter markings
- Fluorescent flexible tip to confirm placement of catheter through vocal chords
- Smooth tip to minimize trauma during intubation
- Sterile, single-use

Specifications

Description	Endotracheal tube introducer
Part Number	BC-01
NSN	TBD
Packed Size	4in x 28.5in (flat) (102mm x 724mm)
Packed weight70oz (20g)
Marking intervals	10cm to 50cm
Shelf life	3 years
Lubricant	Should be used (not included)

PerSys Medical
5310 Elm St.
Houston, TX 77081 USA
1.888.737.7978
info@ps-med.com

www.persysmedical.com

